


MANCHESTER 1824

The University of Manchester

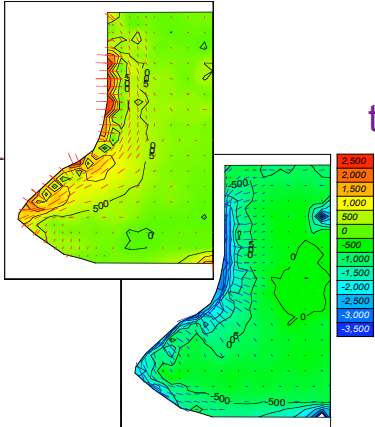

Engineering

Advanced diffraction techniques for Residual Stress determination

Dr Michael Preuss
University of Manchester

Oxford School on Neutron Scattering

1

MANCHESTER 1824

The University of Manchester

Introduction

- Residual stresses in materials
- Principles of measuring residual stresses by diffraction
- Neutron and Synchrotron X-ray diffraction
 - Properties
 - Facilities
- Case Studies / Questions
- From Engineering to Physical Metallurgy – Understanding plasticity
- Conclusions

Oxford School on Neutron Scattering

2

What are residual stresses

Deformation mismatch

Example: Welding




Oxford School on Neutron Scattering


Residual Stresses

- Internal stresses
- Caused by misfit
 - Type I
 - Type II
 - Type III


Bent bar:

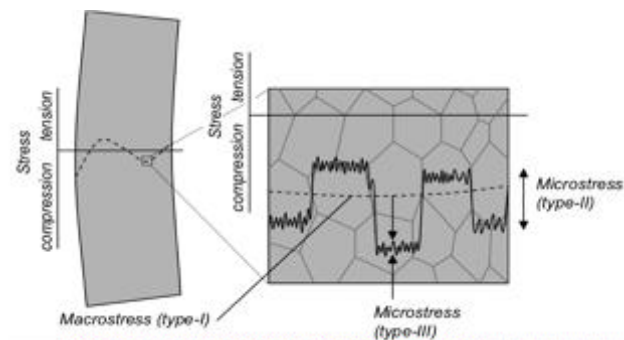

Oxford School on Neutron Scattering

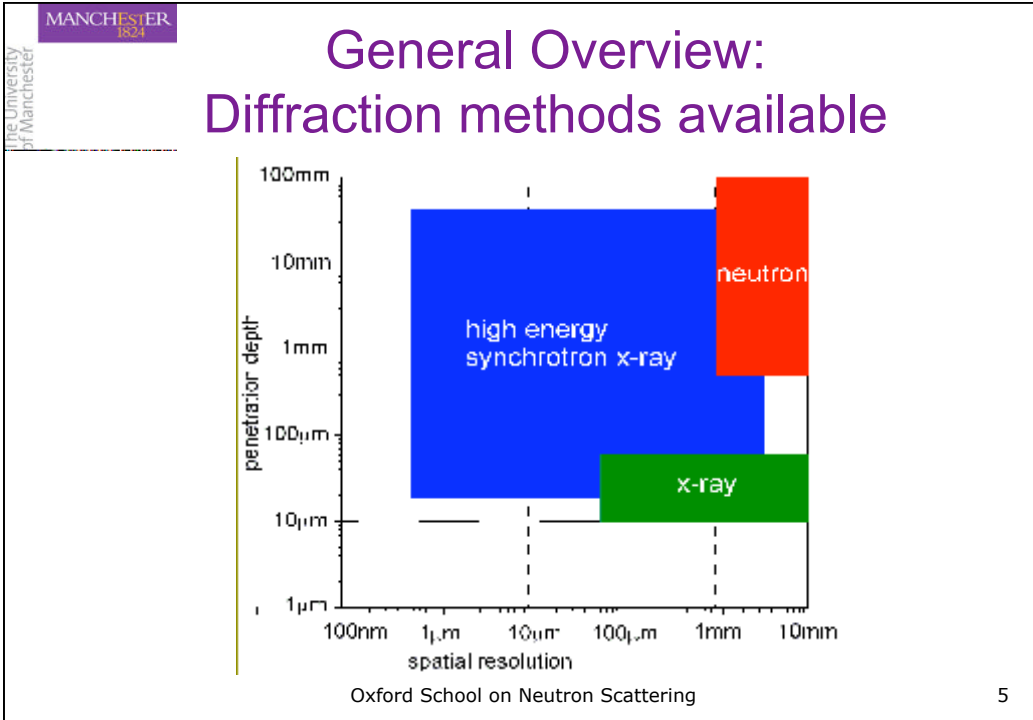

5

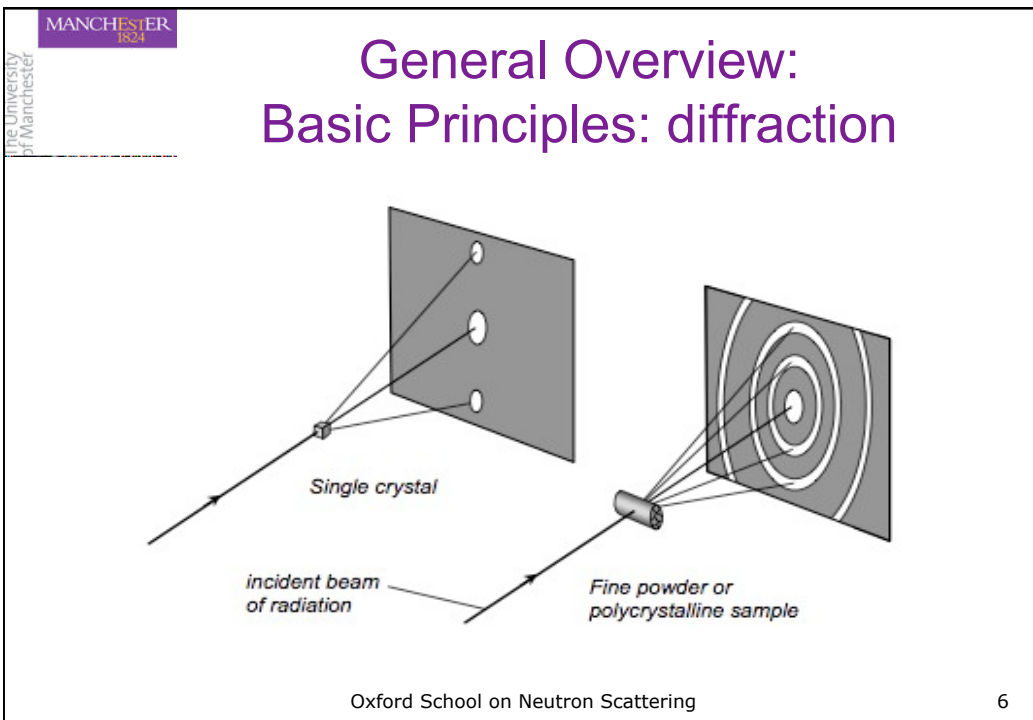

6

MANCHESTER 1824
The University of Manchester

General Overview: Basic Principles: diffraction

$\lambda = 2d \sin \theta$

- Diffraction measures elastic lattice strain as peak shifts
- Uses the poly-crystalline lattice planes as internal strain gauges

Oxford School on Neutron Scattering 7

MANCHESTER 1824
The University of Manchester

Basic Principle

Oxford School on Neutron Scattering 8

MANCHESTER 1824
The University of Manchester

Bragg scattering angle

X-ray wavelength: $\lambda = \frac{hc}{E} = \frac{12.39}{[keV]}$

Bragg's law: $\lambda = 2d \sin \theta$

Scattering Angle

$$\theta = \arcsin\left(\frac{12.39}{2d[keV]}\right)$$

↓

$$\frac{12.39}{[keV]} = 2d \sin \theta$$

Oxford School on Neutron Scattering 9

MANCHESTER 1824
The University of Manchester

d₀ variation

Example of d₀ variation across a tubular Nickel weld

The Vegard Law
Example: Nb in Zr

Oxford School on Neutron Scattering 10

General Overview: Basic Principles

- Measured strains have to be converted into stresses! (Hooke's law)

$$\varepsilon = \frac{a - a_0}{a_0} = \frac{d - d_0}{d_0}$$

- Often requires the unstrained lattice parameter a_0

e.g. isotropic triaxial
along principal
directions:

$$\varepsilon_{11} = \frac{1}{E} [\sigma_{11} - \nu(\sigma_{22} + \sigma_{33})]$$

$$\varepsilon_{22} = \frac{1}{E} [\sigma_{22} - \nu(\sigma_{33} + \sigma_{11})]$$

$$\varepsilon_{33} = \frac{1}{E} [\sigma_{33} - \nu(\sigma_{11} + \sigma_{22})]$$

To calculate a stress direction:

$$\sigma_{11} = \frac{E}{(1 + \nu)(1 - 2\nu)} [(1 - \nu)\varepsilon_{11} + \nu(\varepsilon_{22} + \varepsilon_{33})]$$

(Attention: not always this simple!)


Why do we like neutrons ?

- Part of the nucleus
- Same mass as protons
- Interesting wavelength/mass relationship:

$$\lambda = \frac{h}{p} = \frac{h}{mv}$$

Planck

Mass * Velocity

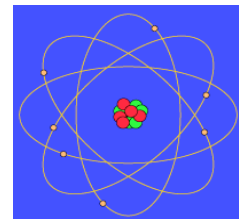

- "Thermal" neutrons: wavelength similar to those of X-rays 0.5-5Å similar to atomic spacing in solids
- Allows cubic gauge volumes!
- Relatively divergent beam !!

MANCHESTER 1824

the University of Manchester

Neutron and Synchrotron Sources


Neutron:

- Reactor Sources (Fission)
 - Constant wavelength/Single Peak
- Accelerator Sources (Spallation)
 - Time-of-flight / Full Spectra / Rietveld


Synchrotron:

- Monochromatic λ and white beam


ESRF
Grenoble, France




ILL
Grenoble, France



ISIS
Chilton, UK



Diamond
Chilton, UK



Oxford School on Neutron Scattering

13


MANCHESTER 1824

the University of Manchester

General Overview: Strain Scanning

Diffracting Gauge Volume: The volume element defined by the incident slits and diffraction slits

Neutron Diffraction

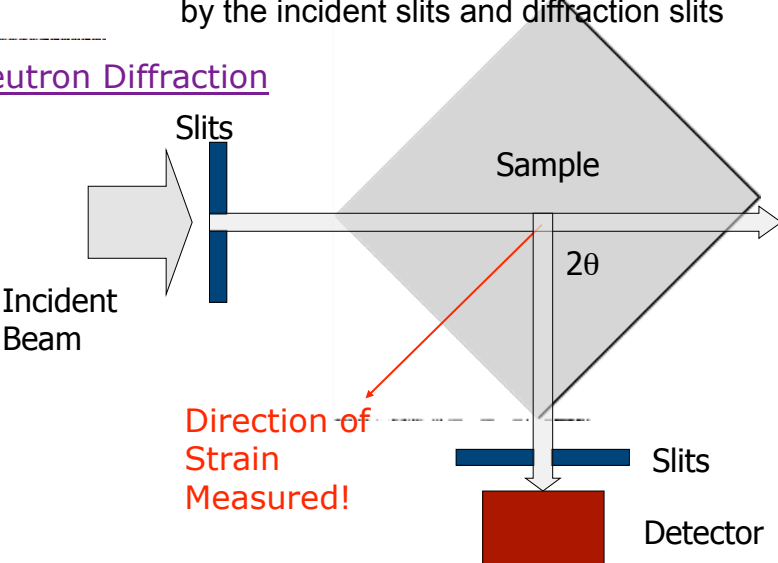

Incident Beam

Slits

Sample

2θ

Direction of Strain Measured!

Slits

Detector

Oxford School on Neutron Scattering

14


In comparison: Synchrotrons

- Synchrotron X-rays are energy X-rays (10-300keV) produced using “synchrotron” accelerators
- Main difference is the wavelength/energy
- Penetration depends on wavelength
- Very high intensity:
 - X-ray tube $\sim 10E8$
 - Synchrotron $\sim 10E15$



In comparison: Synchrotrons

Brilliance of the X-ray beams
(photons / s / mm² / mrad² / 0.1% BW)

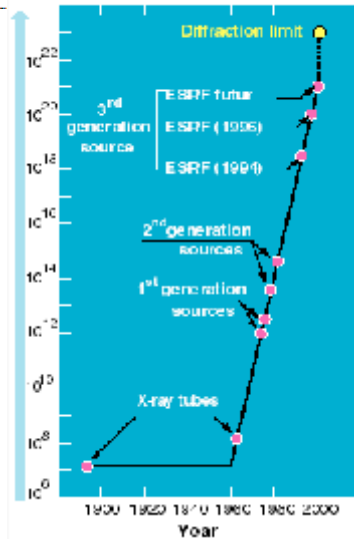

- Synchrotron sources provide very intense (million times more flux than a lab source) high energy beams
- beam is highly parallel (10¹² times more brilliance than a lab source)
- at energies of 40-80keV penetrations of many mm possible
- small micron sized beams

the University of Manchester

MANCHESTER 1824

How does a Synchrotron work ?

- X-rays produced by sending high energy (9GeV) electrons round a ring
 - radial acceleration causes emission of electromagnetic radiation
 - low energy radiation
 - much greater radiation if you insert devices which bend beam sharply
- X-ray beam produced by bending magnets, undulators and wigglers

Oxford School on Neutron Scattering

17

the University of Manchester



MANCHESTER 1824


How does a Synchrotron work ?

- Bending magnets create a wide spectrum of X-ray radiation
- Wigglers are more intense because bend beam many times
- Undulators bend the beam such that radiation interferes to create very high fluxes of certain wavelengths (determined by spacing and number of magnets)
- Highly parallel beam is produced

Oxford School on Neutron Scattering

18

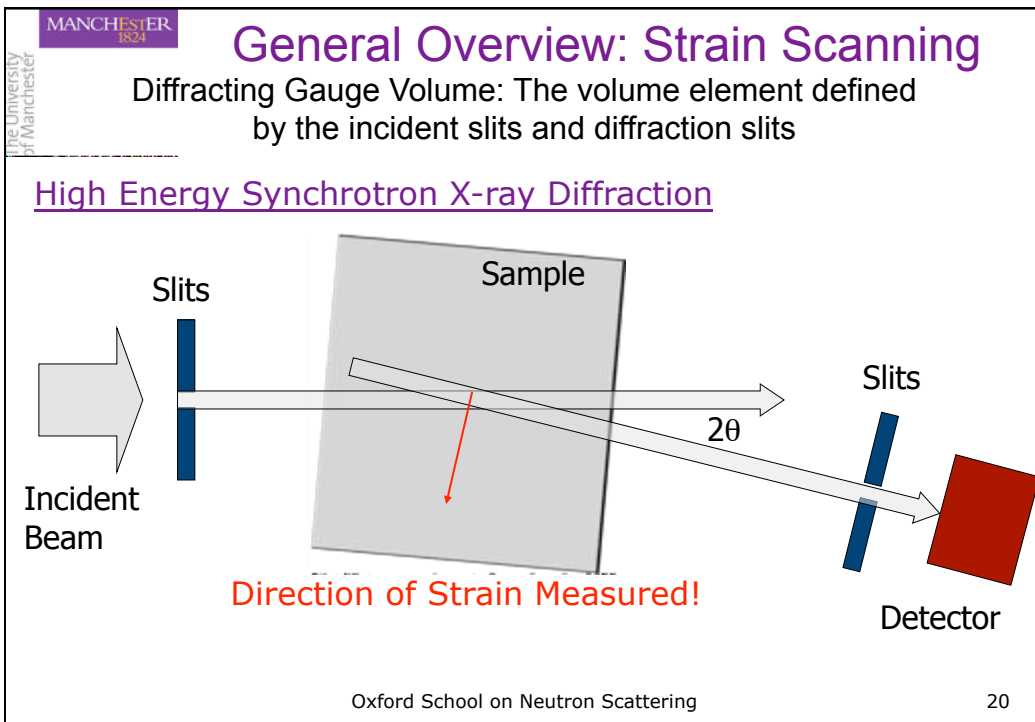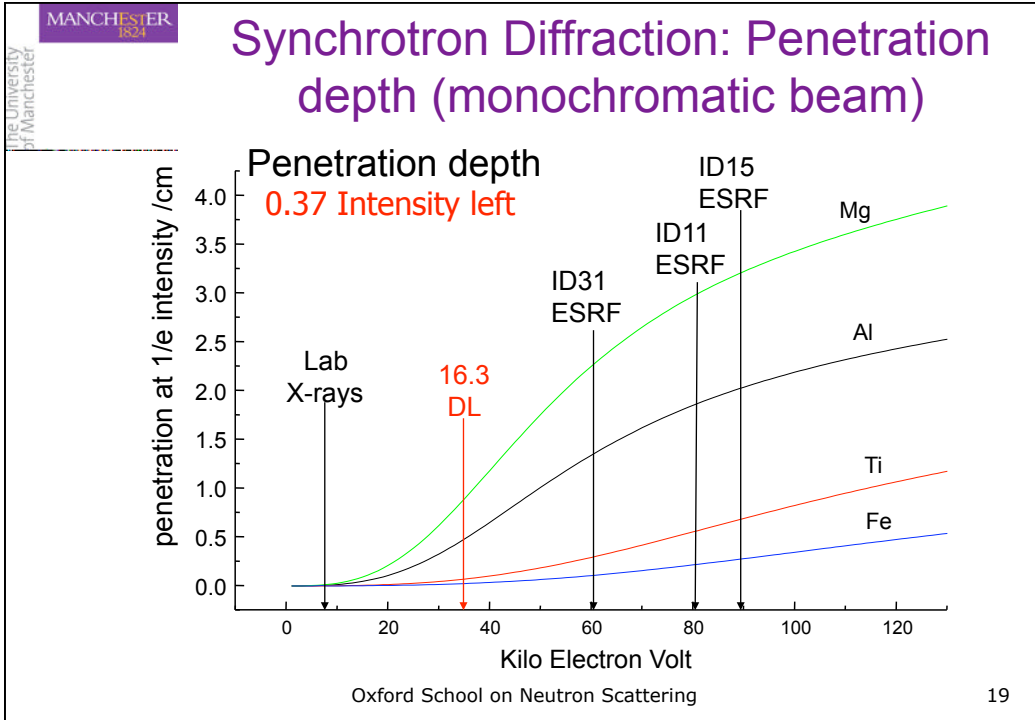

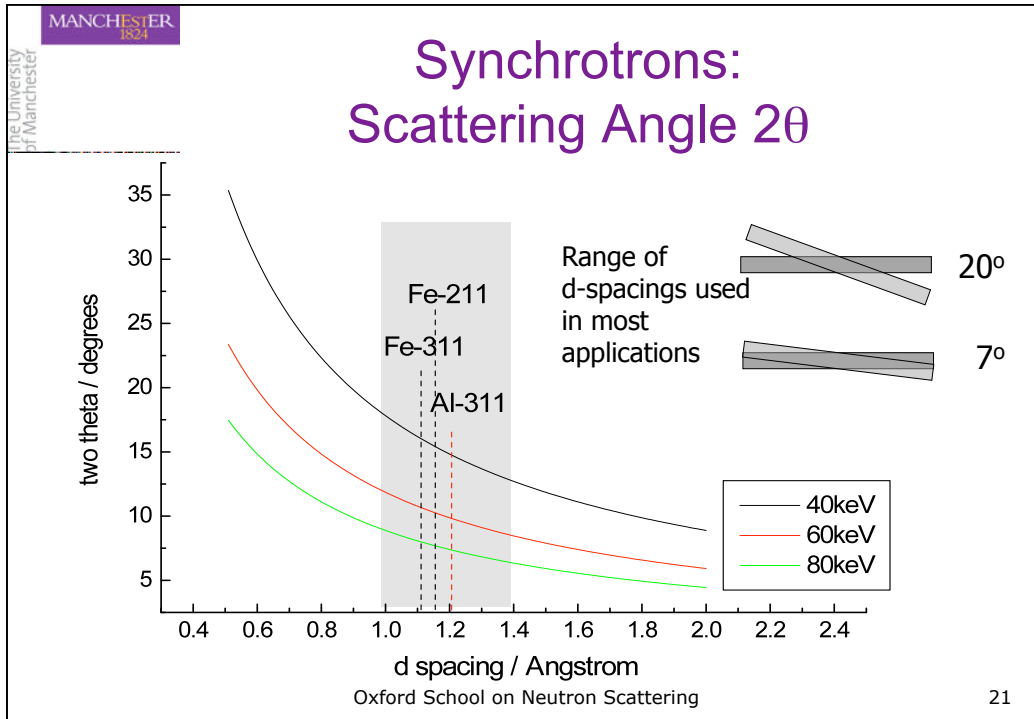

MANCHESTER 1824
the University of Manchester

General Overview: Diffracting Gauge Volume

Volume element of the material in which the recorded scattering takes place

- Results in averaged d-spacing (powder diffraction - many grains)
- Defines the minimum spatial resolution of the method (around 1mm^3 minimum gauge volume when using neutron diffraction)
- and type of residual stress resolved (macro-stress or type-I usually. Type-II for two phase materials).
- Use the largest possible gauge volume for your specific issue in order to minimise counting time

Oxford School on Neutron Scattering

22

Near surface measurements

Neither peak shift (strain) nor measurement location is correct near a surface!

- Partial filling of sampling gauge gives a peak shift - **need to correct peak shift**
- Translator records centre of gauge which is rarely the centre of gravity of diffracting region
 - **need to correct gauge position**

Neutron Properties

- Neutrons are scattered by atomic nuclei (electrons and X-rays which are scattered by the electron cloud).
- Since the scattering is nuclear process, scattering amplitude varies greatly for different isotopes of same element and in a unpredictable manner from element to element. X-ray and electron scattering increase monotonically with atomic number

MANCHESTER 1834

the University of Manchester

Neutron Properties

- Random Scattering length
- Penetration depth independent of energy/wavelength
- Electrically neutral
- Great penetration
- Low flux/intensity

Economic Depth	Al	Steel	Cu	Ti	Ni	SiC
mm	250	37	40	27	24	200

Oxford School on Neutron Scattering 25

MANCHESTER 1834

the University of Manchester

Neutron Scattering

Coherent XS:

$$I_{trans}(t, \lambda) = I_{inc}(\lambda) e^{-\Sigma^* t}$$

Coherent XS ~ Signal

Incoherent XS ~ Background

Absorption XS ~ 1/Intensity

}

Penetration depth ~


1/Sum of Scatt XS

Oxford School on Neutron Scattering 26

MANCHESTER 1834
the University of Manchester

Research Reactors

- Fission in Reactor Core
 - Moderated neutrons
 - Monochromators in guide
- **“Constant Wavelength”**
- Many Facilities in Europe:
 - ILL, SINQ, FRM-2 (G), Petten (NL), ...
 - Generally low flux except ILL and FRM-2




Oxford School on Neutron Scattering 27

MANCHESTER 1834
the University of Manchester

Single Wavelength at Reactor


Single-wavelength instrument: D1A at the ILL
New instrument at ILL: SALSA

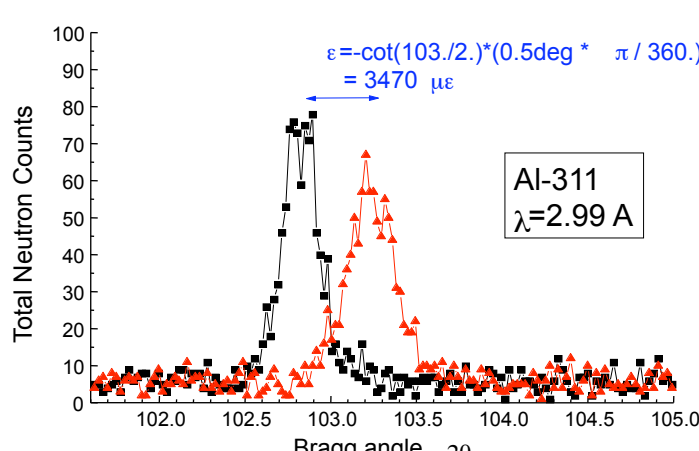

Oxford School on Neutron Scattering 28

Which peak gives us the pure macrostress response ?


In-situ Loading on a neutron diffraction beam line

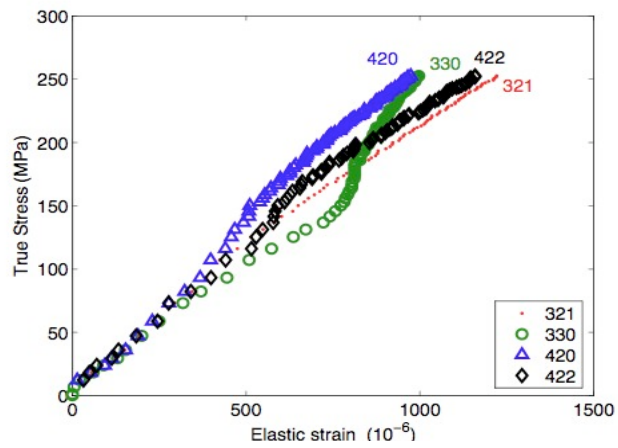

Oxford School on Neutron Scattering

29

SALSA, ILL

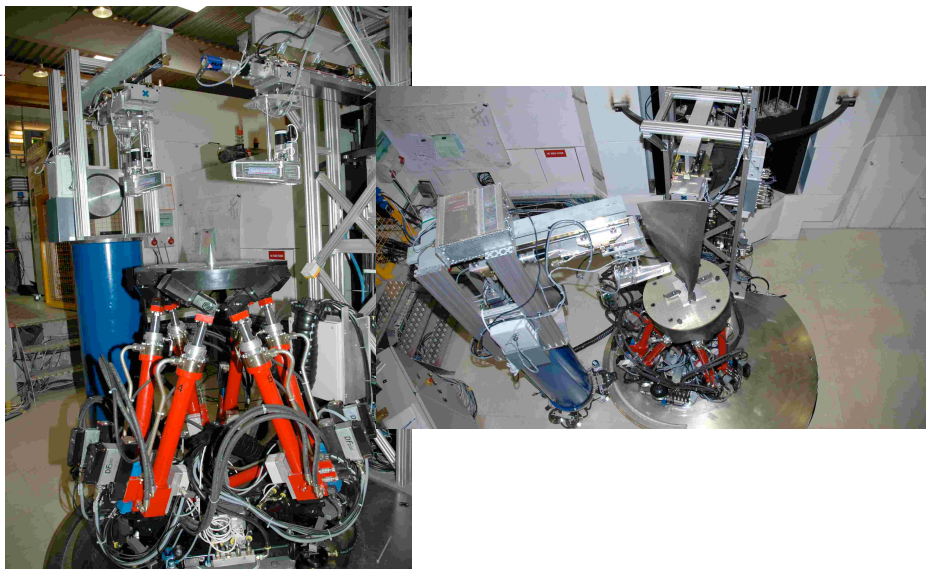

Oxford School on Neutron Scattering

30

MANCHESTER 1824

the University of Manchester

Typical Diffractometer at Synchrotron (here ID31 at the ESRF)

Detector

Sample Stage

Oxford School on Neutron Scattering

31

MANCHESTER 1824

the University of Manchester

Diamond gauge volume

Normal

In-plane

Two orthogonal strain components sample different material volume

Analyscr crystal for partly filled Gauge volume necessary!

Oxford School on Neutron Scattering

32

MANCHESTER 1824
the University of Manchester

Time of flight method

- Sharp pulse leaves source
- High energy neutrons (short λ) travel faster and arrive first, low energy (long λ) last

$$\lambda = ht/ml$$
 where l is the path length and t time of flight
- a single stationary detector records whole diffraction spectrum as a function of time of flight
- neutrons travel at $\sim 100\text{m/s}$ (speed of sound)

$$\lambda = 2d \sin \theta \quad \text{with } \theta \text{ fixed} \quad \text{i.e. } \lambda \text{ proportional to } d$$

Oxford School on Neutron Scattering 33

MANCHESTER 1824
the University of Manchester

Spallation Sources: Time of Flight

Detector

Flight path L

Sample

Fast neutrons arrive earlier at detector!

Neutron Pulse

$$\lambda = \frac{h}{p} = \frac{h}{mv} = \frac{ht}{mL}$$

+

$$\lambda = 2d \sin \theta$$

=

Time-of-Flight:

$$d = \frac{h}{2mL \sin \theta} t$$

Oxford School on Neutron Scattering 34

MANCHESTER 1824
The University of Manchester

Spallation Sources: Measurement of Strain

Strain: $\varepsilon = \frac{a - a_0}{a_0} = \frac{\lambda - \lambda_0}{\lambda_0} = \frac{t - t_0}{t_0}$


Cubic gauge
volume !

↑

Time-of-Flight:

$$\lambda = \frac{h}{2mL \sin \theta} t$$

Fixed




ENGIN-X at ISIS

Oxford School on Neutron Scattering 37


MANCHESTER 1824
The University of Manchester

ENGIN-X, ISIS

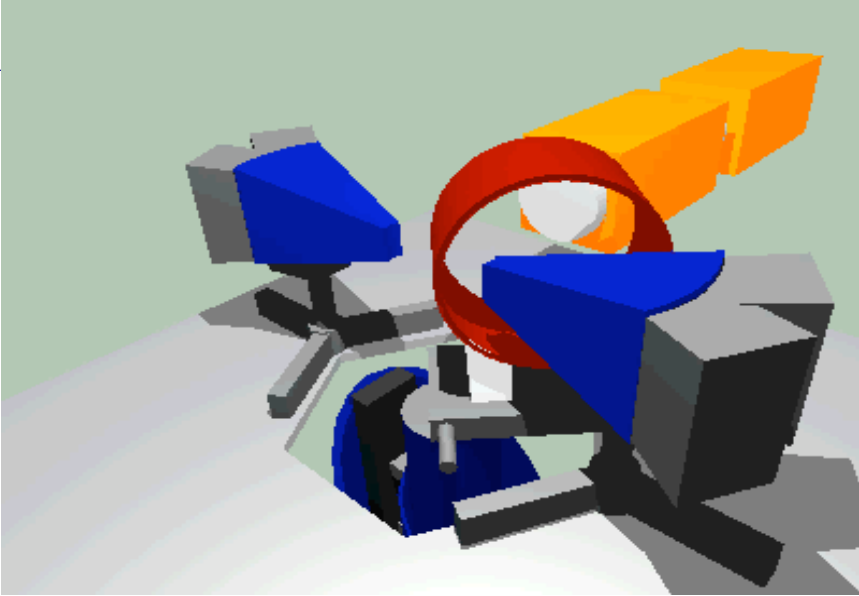

Oxford School on Neutron Scattering 38


ENGIN-X, ISIS

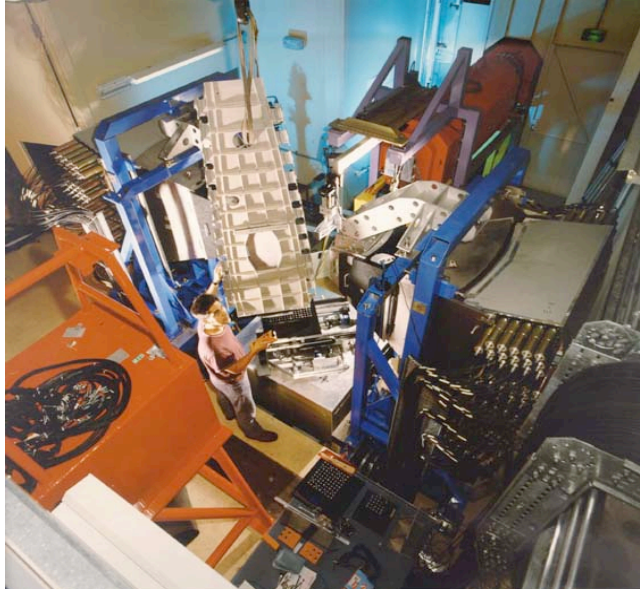

Oxford School on Neutron Scattering

Energy Dispersive Synchrotron Diffraction

- Larger range of wavelengths available
- Energy/Strain Resolution up to $10E-5$
- Higher penetration depth
- More elongated GV

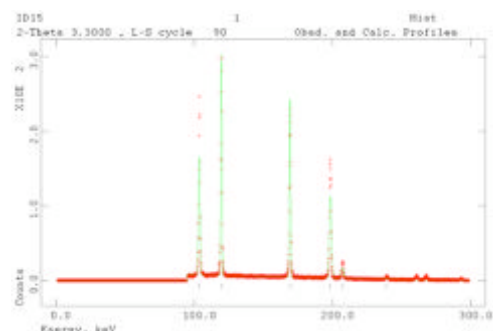

$$\lambda = \frac{hc}{E} = \frac{12.39}{[keV]}$$



$$\epsilon = \frac{E_0}{E} - 1 \quad \text{Strain}$$

MANCHESTER 1824

the University of Manchester

Case Study: Inertia Friction Welding

Solid state joining of compressor, turbine discs and shafts

10th Oxford School on Neutron Scattering

41

MANCHESTER 1824

the University of Manchester

Case Study: Inertia Friction Welding

10th Oxford School on Neutron Scattering

42

MANCHESTER 1824
The University of Manchester

How would you measure such a sample ?

143mm diameter test inertia friction welds




Oxford School on Neutron Scattering

43

MANCHESTER 1824
The University of Manchester

Hoop stresses in IFW'd nickel-base superalloy

All stress units are in MPa

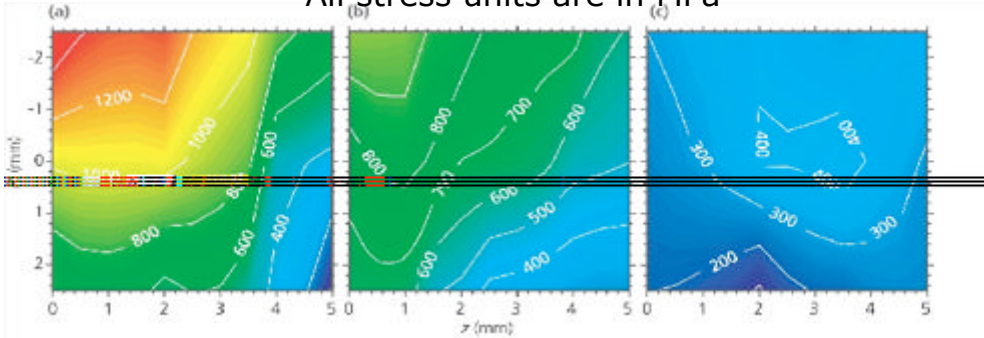

As-welded Conventional PWHT Modified PWHT

Residual stress measurements were used to develop a new PWHT

Oxford School on Neutron Scattering

44

MANCHESTER 1824
The University of Manchester

Railway Rails

- Slices were cut from the rail to measure the horizontal and vertical stresses. Longitudinal stresses were lost
- Measurements were carried out by using neutron and synchrotron x-ray diffraction

Stress (MPa)

-400 -300 -200 -100 0 100 200 300 400

Oxford School on Neutron Scattering 45

MANCHESTER 1824
The University of Manchester

Case Study: Strain Mapping of a TIG weld

2D Map of Residual Strain about the End of a TIG Weld at 100 μ m Resolution


This map include 20,00 measurements and took 8 hours to acquire

8mm wide Tig Weld in Aluminium Plate


Oxford School on Neutron Scattering 46

From Engineering to Physical Metallurgy

Single Crystal Anisotropy

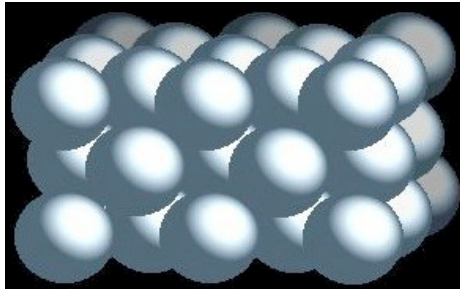

Al, fcc

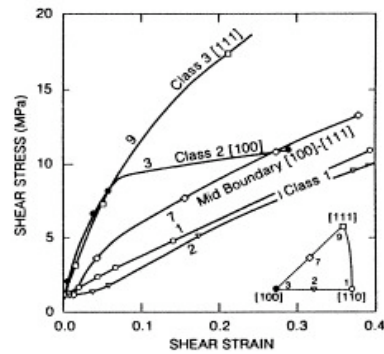


Single Crystal deformation

Deformation heterogeneity


- Polycrystalline deformation is heterogeneous
- Single crystal elastic and plastic anisotropy
- Grain incompatibility during deformation results in intergranular stresses

In-situ loading experiments

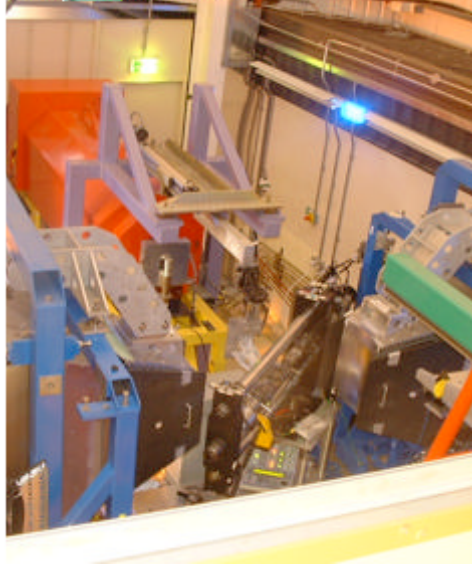


Oxford School on Neutron Scattering

Principle of in-situ loading

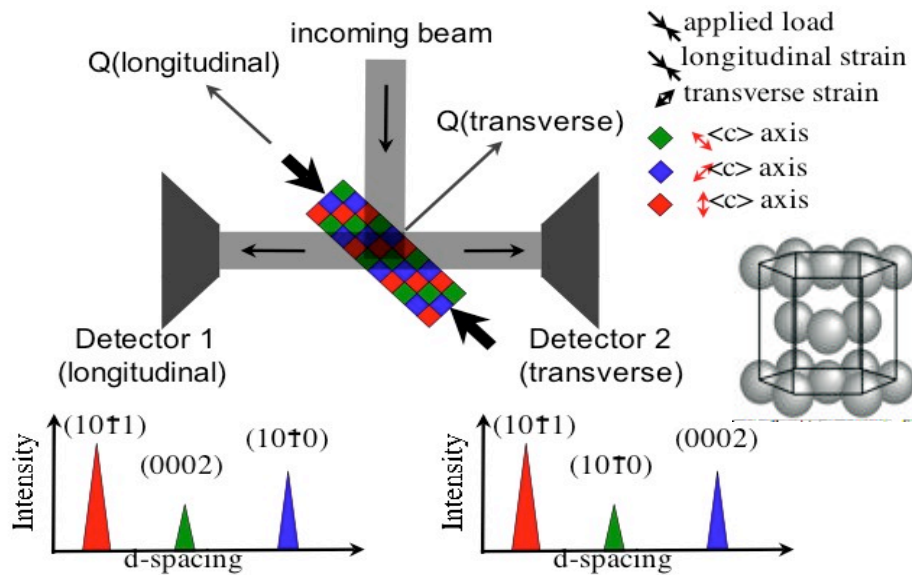


Oxford School on Neutron Scattering

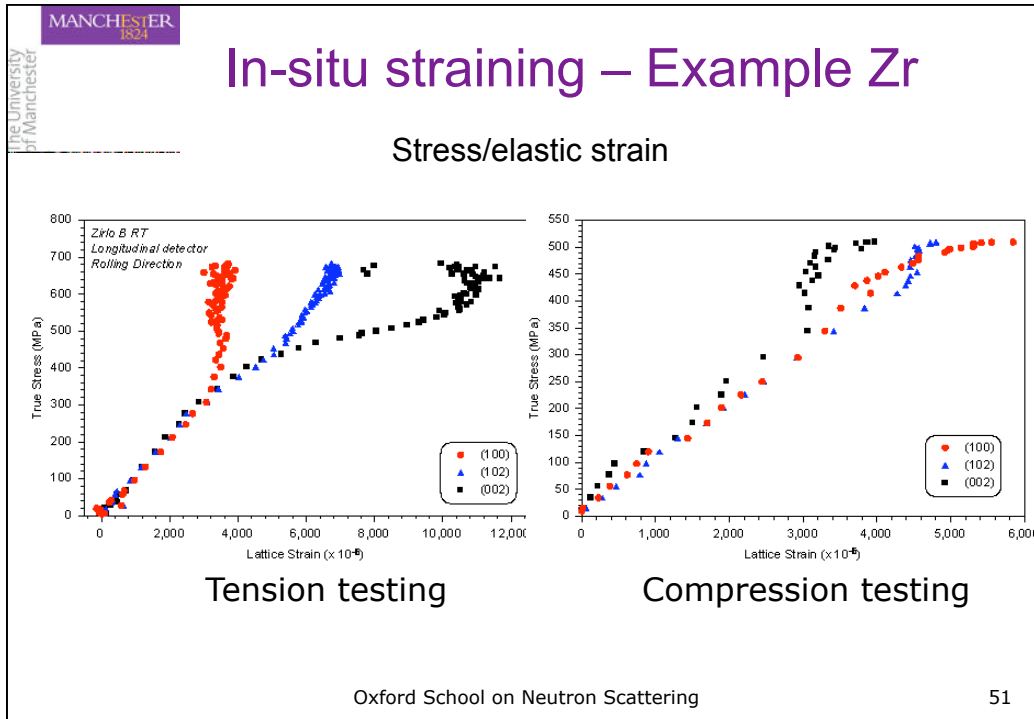

- MANCHESTER 1834
- the University of Manchester
- ## Modelling deformation
- Micromechanics
 - Dislocations, particles, grain boundaries (grain size), interstitial atoms
 - Continuum mechanics:
 - Stresses and strains
 - Intergranular stresses
 - Polycrystal plasticity
 - Mean field methods, i.e. every grain has the same matrix
 - Finite element methods
 - Each grain has a characteristic neighbourhood
 - Predict maximum and minimum stresses ?
- Oxford School on Neutron Scattering 52

EPSC Modelling


- The elasto-plastic self-consistent model (EPSC), is based on the Eshelby-Hill formulation.
- An elliptical inclusion in an infinite medium.
- The surrounding medium is the average of all orientations.
- The inclusion has uniform stress and anisotropic properties i.e. different orientations have different elastic moduli and plastic deformation is only allowed on specified slip planes.
- The model is capable of simulating multiple thermo-mechanical processes.

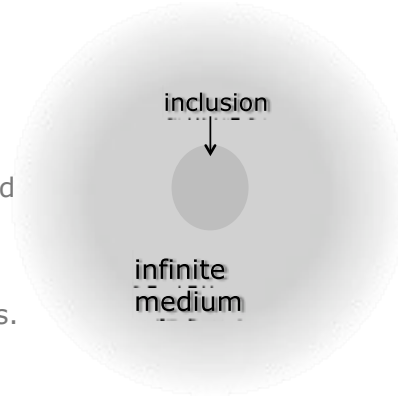

Oxford School on Neutron Scattering


CPFEM

- CPFEM is more computer intensive than EPSC modelling, however, it enables the simulation of specified grain structures.

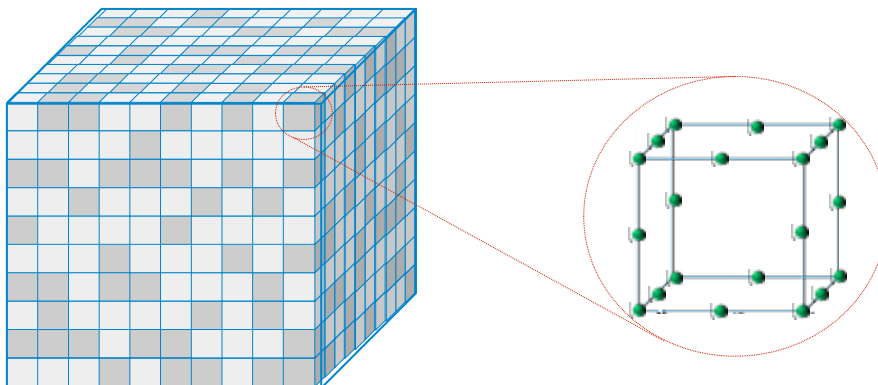

Oxford School on Neutron Scattering

54

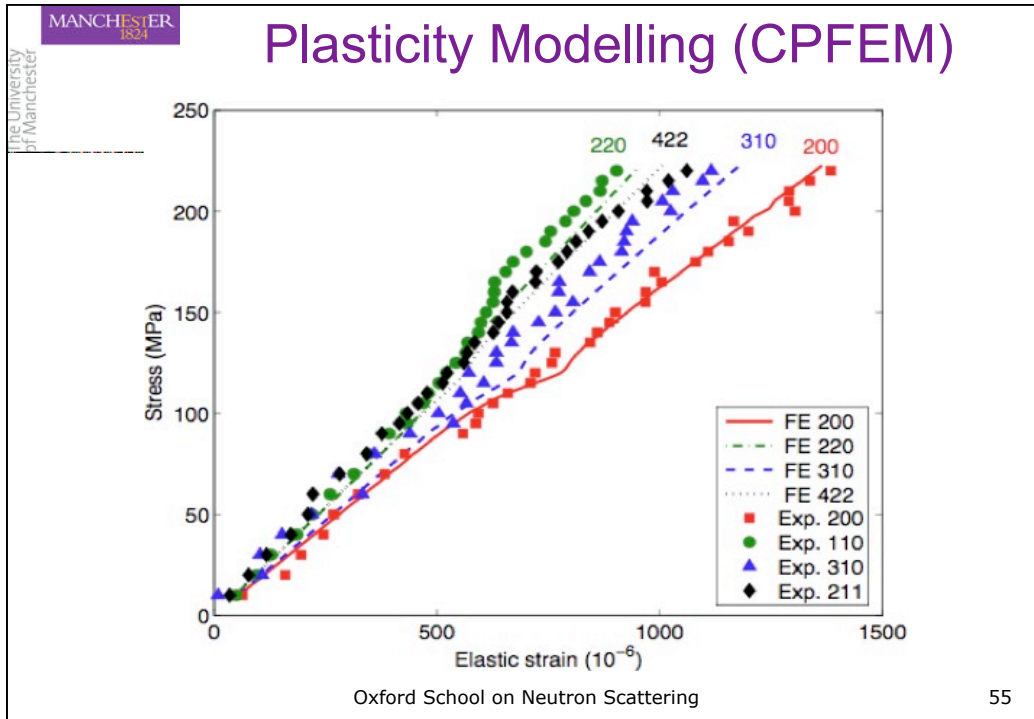

MANCHESTER 1824

the University of Manchester

Attempted General Guidelines: Neutrons

Neutrons:

- Non-destructive, full stress analysis because of cubic Gauge Volume (think three directions)
- Good penetration depth due to neutrality
- Big bulky sample with low stress gradients
- Reasonable spatial resolution independent of atomic number
- Steels, aluminium, nickel, copper zinc or related
- Sample in harsh environment: furnace, cryo. etc.
- Phase analysis with Rietveld analysis

Not-so good: near surface or thin materials, titanium, boron cadmium, fast, high-spatial resolution, high instrumental resolution, hydrogenous materials

Oxford School on Neutron Scattering

56

Attempted General Guidelines: Synchrotrons

Synchrotrons:

- Non-destructive, fast strain mapping, mostly single peak
- Light alloys (small atomic number)
- High spatial resolution aluminium-titanium (think microns)
- High instrumental resolution (small peak width)
- Near surface measurement because of analyser crystal
- Bulk materials / larger atomic number with energy-dispersive method
- Polymers

Not so good at: Steels and higher, big bulky samples, harsh environments, diamond shaped GV